

DISPOSABLE HEALTHCARE PRODUCTS



Surgical Suction

WORKING IN PARTNERSHIP WITH
HEALTHCARE PROFESSIONALS

Pennine[®]
HEALTHCARE



Contents

Suction Connecting Tubing	Suction Connecting Tubing	3
	Soft Suction Connecting Tubing	4
	Conductive Tubing with Anti-static Line	4
	Flexi-Rib Tubing	4
	Tubing with Capped Vacuum Control	4
	Capped Vacuum Control Connector	4
	5-in-1 Connector	4
Link Suction Connecting Tubing	Link Suction Connecting Tubing	5
	Box Link Suction Connecting Tubing	5
Suction Devices	Link Yankauers	6
	Link Yankauer Sets	7
	Derby Yankauers	8
	Selectovac Yankauer	8
	Poole Suckers	8
	Suction Probes	8
Poole Suckers	Derby Yankauer Sets	9
	Poole Sucker Sets	9
	Poole Suction & Flexi-Rib Tubing Sets	9
	Suction Probe Sets	9
Suction Probes	Micro-Flex and Neuro-Flex Suction Probes	10
	Endo-Flex Suction Devices	10
	Sigmo-Flex Suction Devices	10
	Frazier Probes	10
Insufflation	Insufflation Tubing	11
Single-use Medical Devices	Easyflow, Haemo-Ligator & Atiomo DyeSeal	12
Uterine Aspiration Curette System	Aspiration Curettes	13
	Karman Curettes	13
	Wide Bore Aspiration Tubing (Flexi-Cuff)	13
Hy-G-Vac Suction Collection System	Hy-G-Vac Suction Collection System	14
Accessories	Sharps Pad	15
	Tip Cleaner	15
	Diathermy Pencil	15
	Camera Drape	15



Suction Connecting Tubing

A complete system to connect suction sources to suction waste collection systems, suction catheters, yankauers, suction probes and other suction devices. Pennine offer a range of easily connected, interchangeable products in a choice of tubing length and diameter.

Sterile double wrapped individual peel pouches.

'F' (Female) connector fits suction machines and male tapered connections for instrumentation.

'M' (Male) connector push fits into the 'F' connector and is normally used to connect to the standard funnel fitted on suction catheters. Can also be used to connect lengths of tubing together for longer lengths.

Application of 'M' connector with 'VC' (vacuum control) vent converts a standard catheter to vacuum control application.

F Female
FF Female/Female Connectors
FFM Female/Female/Male Connectors
FFM/VC Female/Female/Male Connectors with Vacuum Control



Product Code	Type	Length	CH	Inside Diam. mm	Qty. Per Box
CT-4312	F/F	1.8m	24CH	5	25
CT-4362	F/F	2.7m	24CH	5	20
CT-4372	F/F	3.5m	24CH	5	20
CT-4314	F/F/M	1.8m	24CH	5	25
CT-4364	F/F/M	2.7m	24CH	5	20
CT-4334	F/F/M	3.0m	24CH	5	20
CT-4394/VC	F/F/M/VC	2.7m	24CH	5	20
CT-4374/VC	F/F/M/VC	3.5m	24CH	5	20
CT-4314/VC	F/F/M/VC	1.8m	24CH	5	25
CT-4324/VC	F/F/M/VC	2.0m	24CH	5	25
CT-4334/VC	F/F/M/VC	3.0m	24CH	5	20
CT-4012	F/F	1.8m	26CH	6	25
CT-4022	F/F	2.0m	26CH	6	25
CT-4042	F/F	2.7m	26CH	6	20
CT-4032	F/F	3.0m	26CH	6	20
CT-4072	F/F	3.5m	26CH	6	20
CT-4092	F/F	4.0m	26CH	6	15
CT-4062	F/F	6.0m	26CH	6	10
CT-4020	F/open end	2.0m	26CH	6	25
CT-4030	F/open end	3.0m	26CH	6	20
CT-4014	F/F/M	1.8m	26CH	6	25
CT-4024	F/F/M	2.0m	26CH	6	25
CT-4044	F/F/M	2.7m	26CH	6	20
CT-4034	F/F/M	3.0m	26CH	6	20
CT-4064	F/F/M	6.0m	26CH	6	10
CT-4014/VC	F/F/M/VC	1.8m	26CH	6	25
CT-4024/VC	F/F/M/VC	2.0m	26CH	6	25
CT-4034/VC	F/F/M/VC	3.0m	26CH	6	20
CT-4622	F/F	2.0m	30CH	8	10
CT-4632	F/F	3.0m	30CH	8	10
CT-4672	F/F	3.5m	30CH	8	10
CT-4642	F/F	4.0m	30CH	8	10
CT-4624	F/F/M	2.0m	30CH	8	10
CT-4634	F/F/M	3.0m	30CH	8	10
CT-4674	F/F/M	3.5m	30CH	8	10
CT-4634/VC	F/F/M/VC	3.0m	30CH	8	10



- Softer PVC grade tubing
- Not to be used in high pressure applications

Soft Suction Connecting Tubing

Product Code	Type	Length	CH	Inside Diam. mm	Qty. Per Box
CTS-4332	F/F	3.0m	24CH	5	20
CTS-4022	F/F	2.0m	26CH	6	25
CTS-4032	F/F	3.0m	26CH	6	20
CTS-4072	F/F	3.5m	26CH	6	20
CTS-4092	F/F	4.0m	26CH	6	15
CTS-4640	F/open end	4.0m	30CH	8	10

- Anti-static line, prevents discharge



Conductive Tubing with Anti-static Line

CCT-4112	F/F	1.8m	26CH	6	25
CCT-4122	F/F	2.0m	26CH	6	10
CCT-4142	F/F	2.7m	26CH	6	25
CCT-4132	F/F	3.0m	26CH	6	20
CCT-4172	F/F	3.5m	26CH	6	20

- Ribbed connecting tubing for enhanced drapability
- Non-slip outer rib reduces friction
- Resistant to coiling, kinking
- Maximum pressure 500mm HG Mercury
- Colour coded funnel identification
- Multi-fit to all suction devices and sources
- Economical box size



Flexi-Rib Tubing

CT-4922	F/F	2.0m	26CH	6	60
CT-4932	F/F	3.0m	26CH	6	50
CT-4972	F/F	3.5m	26CH	6	50
CT-4924	F/F/M	2.0m	26CH	6	60
CT-4934	F/F/M	3.0m	26CH	6	50
CT-4722	F/F	2.0m	28CH	7	60
CT-4732	F/F	3.0m	28CH	7	50
CT-4772	F/F	3.5m	28CH	7	40
CT-4724	F/F/M	2.0m	28CH	7	60
CT-4734	F/F/M	3.0m	28CH	7	50
CT-4774	F/F/M	3.5m	28CH	7	40
CT-4714/VC	F/F/M/VC	1.8m	28CH	7	25
CT-4741/VC	F/F/M/VC	2.5m	28CH	7	20

Tubing with Capped Vacuum Control

CT-4015	F/VCU2	1.8m	26CH	6	25
CT-4025	F/VCU2	2.0m	26CH	6	25

Capped Vacuum Control Connector

VCU2	Sterile	100
VC2/Bulk	Bulk, non Sterile	1200

- Available sterile and non-sterile or in bulk quantities



- Clear PVC, single use only
- Can be used to extend length of tubing i.e. to connect two tubes together of either the same or different diameters
- Available either bulk, non-sterile or packed sterile
- Designed for use with all types of tubing

5-in-1 Connector

DP-501	Sterile	100
DP-501/Bulk	Bulk, non sterile	500




This unique system will guarantee an easy and faultless connection between all the individual products in the Link range. Link Tubing is also referred to as Bubble Tubing.

Although the Link products are optimised to be used within the Link system, they can also be easily connected to other products in our range.

- Sterile, double wrapped transparent individual peel pouches
- Disposable PVC tubing
- Available in two diameters, 5mm & 7mm

Available either in plain PVC or with an anti-static black line to prevent static discharge (also known as conductive tubing)

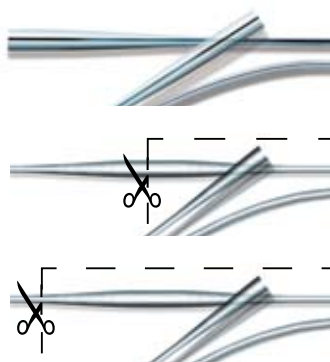

For integral funnel connector - cut here

For integral tapered connector - cut here

Link Suction Connecting Tubing

The Pennine Link system is used for connecting suction sources to suction waste collection systems, suction catheters, yankauer suckers, suction probes and other suction devices. Tubing extruded with an integral female funnel and a tapered male end to provide adaptable connections. The tube can be modified by cutting with scissors.


Product Code	Type	Length	CH	Inside Diam. mm	Qty. Per Box
CT-4513	Plain F/F	1.8m	24CH	5mm	25
CT-4543	Plain F/F	2.7m	24CH	5mm	20
CT-4511	Plain M/F	1.8m	24CH	5mm	25
CT-4541	Plain M/F	2.7m	24CH	5mm	20
CT-4713	Plain F/F	1.8m	26CH	7mm	25
CT-4743	Plain F/F	2.7m	26CH	7mm	20
CT-4711	Plain M/F	1.8m	26CH	7mm	25
CCT-4513	Conductive F/F	1.8m	24CH	5mm	25
CCT-4543	Conductive F/F	2.7m	24CH	5mm	20
CCT-4583	Conductive F/F	3.7m	24CH	5mm	15
CCT-4512	Conductive M/F	1.8m	24CH	5mm	25
CCT-4542	Conductive M/F	2.7m	24CH	5mm	20
CCT-4582	Conductive M/F	3.7m	24CH	5mm	15
CCT-4713	Conductive F/F	1.8m	26CH	7mm	25
CCT-4743	Conductive F/F	2.7m	26CH	7mm	20
CCT-4783	Conductive F/F	3.7m	26CH	7mm	10
CCT-4712	Conductive M/F	1.8m	26CH	7mm	25
CCT-4742	Conductive M/F	2.7m	26CH	7mm	20
CCT-4782	Conductive M/F	3.7m	26CH	7mm	10

Box Link Suction Connecting Tubing

LINK tubing available in 30 and 50 metre lengths supplied in convenient dispenser boxes. Tube extruded with a bubble at every 0.90m intervals. Simply cut through the bubble at the required length to create either a male or female end which will fit a suction tip or suction catheter.

Product Code	Type	Length	CH	Inside Diam. mm	Qty. Per Box
SL-8005	Plain	30m	24CH	5mm	1 Roll
SL-8055	Plain	50m	24CH	5mm	1 Roll
SL-8056	Plain	50m	26CH	6mm	1 Roll
SL-8007	Plain	30m	28CH	7mm	1 Roll
SL-8205	Conductive	30m	24CH	5mm	1 Roll
SL-8255	Conductive	50m	24CH	5mm	1 Roll
SL-8256	Conductive	50m	26CH	6mm	1 Roll
SL-8207	Conductive	30m	28CH	7mm	1 Roll

- LINK tubing available in 30 and 50 metre lengths
- Three diameters available: 5mm, 6mm and 7mm
- Supplied in convenient dispenser boxes

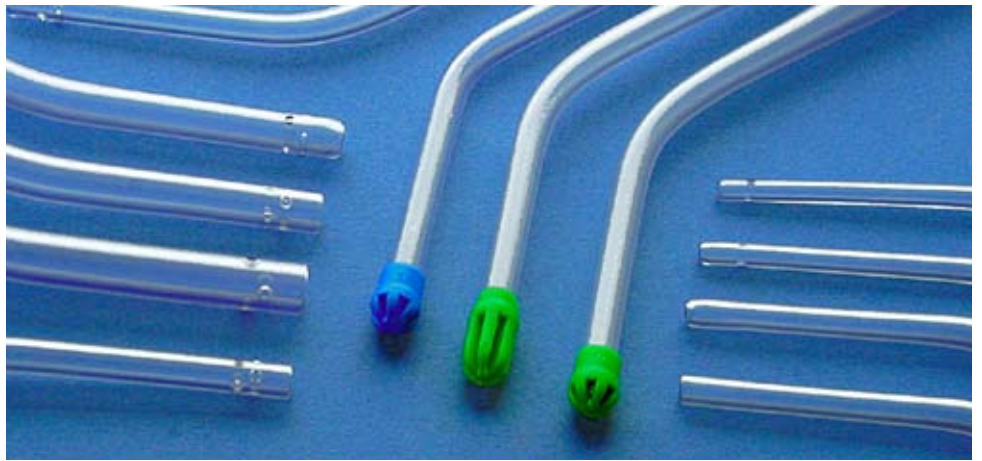


Link Yankauers

Rigid, transparent suction devices for use with suction connecting tubing. For general purpose surgical suctioning. Also available as a set, pre-attached to 2 or 3 metre tubing.

- Sterile, double wrapped transparent individual peel pouches
- Available with or without vacuum control (VC) vent
- Short mini, Mini fine, Midi and Maxi Yankauers available with smooth rounded end or **Super Tip**
- Midi Yankauer available with additional **Atrauma Tip** to reduce tissue grab

Product Code	Type	Use	Length	CH	Qty. Per Box
LY-3400	Short Mini	Paediatric/dental	190mm	15ch	100
LY-3401	Short Mini VC		190mm	15ch	100
LY-4400	Short Mini Supertip		195mm	15ch	100
LY-4401	Short Mini Supertip VC		195mm	15ch	100
LY-3500	Mini Fine Tip	Fine suctioning (e.g. ENT)	270mm	12ch	50
LY-3501	Mini Fine Tip VC		270mm	12ch	50
LY-4500	Mini Fine Supertip		275mm	12ch	50
LY-4501	Mini Fine Supertip VC		275mm	12ch	50
LY-3600	Midi 30° angle	For general and pharyngeal suction	280mm	18ch	50
LY-3601	Midi VC		280mm	18ch	50
LY-3600/PR	Midi 10° angle		270mm	18ch	50
LY-4600	Midi Supertip		280mm	18ch	50
LY-4601	Midi Supertip VC		280mm	18ch	50
LY-3700	Maxi		For large particles and significant volumes of fluid	275mm	22ch
LY-3701	Maxi VC	275mm		22ch	50
LY-4700	Maxi Supertip	275mm		22ch	50
LY-4701	Maxi Supertip VC	275mm		22ch	50
LY-3800	Maxi Plus	For orthopaedic use	280mm	28ch	50
LY-3801	Maxi Plus VC		280mm	28ch	50
LY-3602	Atrauma tip, short soft	For general suction with minimal tissue grab	280mm	18ch	50
LY-3603	As above with VC		280mm	18ch	50
LY-3606	Atrauma tip, short rigid		280mm	18ch	50
LY-3607	As above with VC		280mm	18ch	50
LY-3608	Atrauma tip, long rigid		280mm	18ch	50
LY-3609	As above with VC		280mm	18ch	50

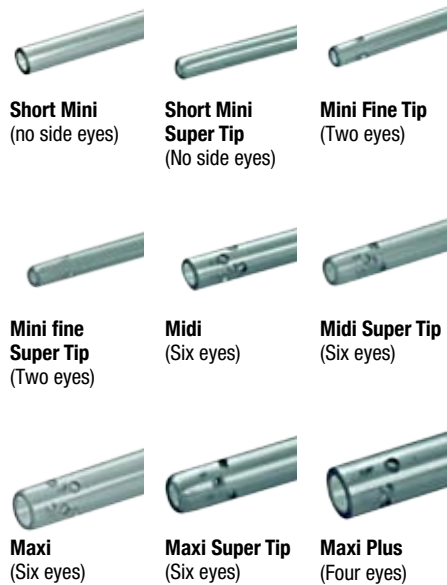

Atrauma Tips





Link Yankauer Sets

Link yankauers pre-connected to Pennine suction connection tubing in 2.0 and 3.0 metre lengths. For general purpose surgical suctioning.

Sets with 2m tubing

- Sterile, double wrapped transparent individual peel pouches
- Available with or without VC (vacuum control) vent

Product Code	Tubing Diam./Length	Qty. Per Box	Yankauer Type
LYS-5420	6mm / 2m	10	LY-3400
LYS-5421	6mm / 2m	10	LY-3401
LYS-5520	6mm / 2m	10	LY-3500
LYS-5521	6mm / 2m	10	LY-3501
LYS-5620	6mm / 2m	10	LY-3600
LYS-5621	6mm / 2m	10	LY-3601
LYS-5720	6mm / 2m	10	LY-3700
LYS-5721	6mm / 2m	10	LY-3701
LYS-6820	7mm / 2m	10	LY-3800
LYS-6821	7mm / 2m	10	LY-3801
LYS-3622	6mm / 2m	10	LY-3602
LYS-3623	6mm / 2m	10	LY-3603
LYS-3626	6mm / 2m	10	LY-3606
LYS-3627	6mm / 2m	10	LY-3607
LYS-3628	6mm / 2m	10	LY-3608
LYS-3629	6mm / 2m	10	LY-3609

Sets with 3m tubing

Product Code	Tubing Diam./Length	Qty. Per Box	Yankauer Type
LYS-5430	6mm/3m	10	LY-3400
LYS-5431	6mm/3m	10	LY-3401
LYS-5530	6mm/3m	10	LY-3500
LYS-5531	6mm/3m	10	LY-3501
LYS-5630	6mm/3m	10	LY-3600
LYS-5631	6mm/3m	10	LY-3601
LYS-5730	6mm/3m	10	LY-3700
LYS-5731	6mm/3m	10	LY-3701
LYS-6830	7mm/3m	10	LY-3800
LYS-6831	7mm/3m	10	LY-3801
LYS-3632	6mm/3m	10	LY-3602
LYS-3633	6mm/3m	10	LY-3603
LYS-3636	6mm/3m	10	LY-3606
LYS-3637	6mm/3m	10	LY-3607
LYS-3638	6mm/3m	10	LY-3608
LYS-3639	6mm/3m	10	LY-3609



Derby Yankauers

Clear, acrylic suction devices, similar to traditional yankauer suckers. With or without sump tip to reduce potential tissue suction during use.

- Sterile, double wrapped in individual peel pouches
- Available with or without VC (vacuum control) vent
- Available with or without sump tip to reduce the risk of tissue damage



Product Code	Type	Length	CH/ Tip	Inside Diam. mm	Qty. Per Box
YS-3004	with sump tip	265mm	24CH	8mm	50
YS-3005	with sump tip and VC	265mm	24CH	8mm	50
YS-3014	straight tip	265mm	24CH	8mm	50
YS-3015	straight tip and VC	265mm	24CH	8mm	50
YS-3000	Midi Yankauer with straight tip	265mm	20CH	6.66mm	50

Selectovac Yankauer

Surgical suction cannula available with bulb tip.

- Switch system allows the operator to manually open and close suction
- Easy connection to Pennine range of tubing



Product Code	Type	Length	CH/ Tip	Inside Diam. mm	Qty. Per Box
YS-3009	On/Off valve	265mm	24CH	8mm	50

Pooler Suckers

Flexible or Rigid sump sucker with removable inner suction probe. Specifically for abdominal surgery, reducing trauma to surrounding tissues.

- Available with atraumatic end
- Available with angled tip

- Removable inner probe



Product Code	Type	Length	CH/ Tip	Inside Diam. mm	Qty. Per Box
YS-3006	Flexible straight tip	300mm	20CH	6.66mm	50
YS-3008	Flexible (angled tip)	300mm	20CH	6.66mm	50

YS-3007	Rigid (White)	300mm	15CH	6.5mm	50
---------	---------------	-------	------	-------	----

Suction Probes

Pennine suction devices manufactured from rigid transparent PVC, for general suction purposes, orthopaedic suctioning or dental work. Easily attached to Pennine suction tubing.

- Standard probes 20ch
- Available in a wider bore, 30ch for orthopaedic use and for use with wide bore suction connection tubing



Product Code	Type	Length	CH/Tip	Qty. Per Box
SP-3800	Suction probe with male connector.	245mm	20CH	50
SP-3801	As above with VC (vacuum control)	245mm	20CH	50
SP-3900	Flat tip probe with male connector	250mm	20CH	50
SP-3901	As above with VC (vacuum control)	250mm	20CH	50
SP-4200	Maxi probe	280mm	30CH	50
SP-4201	As above with VC (vacuum control)	280mm	30CH	50
SP-4300	Flat tip Maxi probe	255mm	30CH	50
SP-4301	As above with VC (vacuum control)	255mm	30CH	50
MP-3800	Dental suction device	120mm	20CH	50

Please note: All angles and lengths are nominal

Tel: 44 (0) 1332 794880 Email: marketing@penninehealthcare.co.uk www.penninehealthcare.co.uk

LATEX
FREE



- All Pennine suction sets come sterile, double wrapped in individual peel pouches
- Available with or without VC (vacuum control) hole

Derby Yankauer Sets

Derby Yankauer Suckers pre-attached to 2m or 3m connecting tube. Available with or without 'sump' tip to reduce the risk of tissue damage during use.

Set with 2m tubing	CH/Tip	Qty	Set with 3m tubing	Qty Per Box	Tube Diam. mm	Yankauer Type
SYS-5820	24CH	20	SYS-5830	10	6mm	YS-3004
SYS-5821	24CH	20	SYS-5831	10	6mm	YS-3005
SYS-5520	24CH	20	SYS-5530	15	6mm	YS-3014
SYS-5521	24CH	20	SYS-5531	15	6mm	YS-3015

Poole Sucker Sets

Flexible 'sump' sucker pre-attached to 2m or 3m connecting tube. For deep abdominal suctioning or where a 'sump' suction action is required.

- Available with angled tip

Set with 2m tubing	CH/Tip	Qty Per Box	Set with 3m tubing	Qty Per Box	Tube Diam. mm	Yankauer Type
SYS-5920	20CH	10	SYS-5930	10	6mm	YS-3006
SYS-5923	20CH	15	SYS-5933	10	6mm	YS-3008
SYS-5720	15CH	10	SYS-5730	10	6mm	YS-3007

Poole Suction & Flexi-Rib Tubing Sets

3m Flexi-Rib tubing pre-attached to Poole suction probe.

- Removable inner probe

Set with 3m tubing	CH/Tip	Qty Per Box	Tube Diam. mm	Poole Type
SYS-5931	20CH	50	6mm Flexi-Rib	YS3006
SYS-5932	20CH	50	6mm Flexi-Rib	YS3008

Suction Probe Sets

Suction probes available as sets pre-attached to 2m or 3m tubing. For general use (by attachment to 6mm diameter standard connecting tubing) or orthopaedic use (by attachment to 30ch, 8mm wide bore tubing).

Set with 2m tubing	CH/Tip	Box Qty	Set with 3m tubing	Qty Per Box	Tube Diam. mm	Yankauer Type
SPS-5020	20CH	25	SPS-5030	20	6mm	SP-3800
SPS-5021	20CH	25	SPS-5031	20	6mm	SP-3801
SPS-5320	20CH	25	SPS-5330	20	6mm	SP-3900
SPS-5321	20CH	25	SPS-5331	20	6mm	SP-3901
SPS-6020	30CH	10	SPS-6030	10	8mm	SP-4200
SPS-6021	30CH	10	SPS-6031	10	8mm	SP-4201
SPS-6320	30CH	10	SPS-6330	10	8mm	SP-4300
SPS-6321	30CH	10	SPS-6331	10	8mm	SP-4301

SP-3800



MP-3800



SP-3900



SP-4200



SP-4300





Micro-Flex and Neuro-Flex Suction Probes

Special PVC disposable probes with efficient prestrol vacuum control device. Designed for neuro surgery, ENT surgery and in applications where controlled suction is required.


Product Code	Use	Length	CH/ Tip	Qty. Per Box
YS-3148/14	For neuro surgery and other intricate suction procedures.	16cm	14ch Ø 4.66mm	50
YS-3128		15cm	12ch Ø 4.00mm	50
YS-3108	For neuro surgery and other intricate suction procedures.	15cm	10ch Ø 3.33mm	50
YS-3088		15cm	8ch Ø 2.66mm	50
YS-3088(Z)	For otological procedures and other intricate surgical applications	11cm	8ch Ø 2.66mm	50
YS-3068		11cm	6ch Ø 2.00mm	50
YS-3058		11cm	5ch Ø 1.66mm	50

- Sterile, double wrapped individual peel pouches.
- Neuro-Flex (Frazier)
- Micro-Flex (Frazier)
- Micro-Flex (Zoellner)

Endo-Flex Suction Devices

Special probes for laryngoscopy, pharyngoscopy and bronchoscopy procedures. All 12ch tubes with varying probe length. Prestrol vacuum control vent with small connector to fit Pennine connection tubing.

Product Code	Use	Length	CH/ Tip	Qty. Per Box
YS-3128/20	Laryngoscopy suction probe	20cm	12ch Ø 4.00mm	10
YS-3128/26	Pharyngeal suction probe	26cm	12ch Ø 4.00mm	10
YS-3128/40	Bronchoscopy suction probe	40cm	12ch Ø 4.00mm	10
YS-3128/60	Bronchoscopy suction probe	60cm	12ch Ø 4.00mm	10

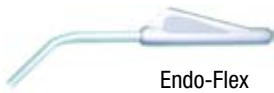

Sigmo-Flex Suction Devices

An Endoflex but with standard non 'VC' vacuum control male connector. 20ch tube for Sigmoidoscopy/Oesophagoscopy.

Product Code	Use	Length	CH/Tip	Qty. Per Box
SP-5800	Sigmoidoscopy suction probe	60cm	20ch Ø 6.66mm	10

Frazier Probes

A single-use probe manufactured from anodized aluminium for extra flexibility. With 'VC' vacuum control hole for precise control over the amount of suctioning. Used in open surgical procedures where delicate suctioning is required.

Product Code	Use	Length	CH/ Tip	Qty. Per Box
MF-1108	For neuro surgery and other intricate suction procedures	230mm	8ch	50
MF-1110		230mm	10ch	50
MF-1112		230mm	12ch	50

- Malleable shaft can be bent to desired position
- Longer probe
- Longer, non-slip handle





Protects the patient from gas borne contamination.

- 0.2µm hydrophobic filters prevents cross contamination from patient to machine and vice versa
- 6mm diameter lightweight tubing in 2.0, 2.5 and 3.0 metre lengths
- Choice of machine end connectors for compatibility with standard insufflators
- Sterile, double-wrapped individual peel pouches

- 1) Open packaging and place in sterile field
- 2) Attach the shorter piece of tubing to the insufflator luer or male barb connection
- 3) Attach the other end to Verres needle or Insufflation port

*INS-1250 For export only

Insufflation Tubing

For use in laparoscopic surgery (minimally invasive) for delivery of clean filtered CO₂. This system is to be used only for the insufflation of carbon dioxide into the abdomen during laparoscopy.

Product Code	Type	Length	Inside Diam. mm	CH/Tip	Qty. Per Box
INS-1320/LL	Tubing with filter, male luer lock/male luer lock	2.0	6.0mm	26CH	10
INS-1330/LL	Tubing with filter, male luer lock/male luer lock	3.0	6.0mm	26CH	10
INS-1330/E	Tubing with filter, open end/funnel connector	3.0	6.0mm	26CH	10
INS-1332	Tubing with filter, male luer lock/push fit connector	3.0	6.0mm	26CH	10
INS-1250	Tubing with filter, male luer lock/Wolf connector*	2.5	6.0mm	26CH	10



INS-1330/LL
INS-1320/LL



INS-1250*



INS-1330/E



INS-1332



Easyflow

Filtered suction system for orthopaedic surgery.

- Suitable for removal of large bone fragments during orthopaedic surgery
- Suitable for hip arthroplasty
- Eliminates blockages
- Suitable for perioperative blood collection
- Rapid filter exchange whilst maintaining suction
- Reduces risk of tissue damage and post operative infections and complications
- Ergonomic handle design ensures comfort and control
- Reversible suction tip provides extra length for intra-femoral suction

Product Code	Use	Contents	Qty. Per Box
EFS-1111	Orthopaedic Surgery	7mm/3m Flexi-Rib Tubing • 2 curved probes • 1 straight probe • Spare filter	10



Haemo-Ligator

Preloaded single-use multi-band suction system. The Pennine Haemo-Ligator is a disposable haemorrhoid suction ligator, ready to use, pre-loaded with five bands and pre-attached to flexible suction tubing.

- Packed clean
- Ergonomic design
- Suitable for use with all commonly available ano-proctoscopes
- Single practitioner procedure

Product Code	Type	Tubing Diam. mm	Qty. Per Box
HLS-001/LF	Haemo-Ligator	7mm / 3m	10



Atiomo DyeSeal

The Atiomo Dyeseal Uterine Manipulator is suitable for all procedures where a Spackman Cannula is used and it has added benefits in Lap & Dye procedures.

- Sterile, single-use
- Cost effective
- No loose parts
- No reprocessing costs
- Reduces risk of cross contamination

Product Code	Type	Qty. Per Box
CW-1200	Uterine Manipulator	10





Aspiration Currettes

Clear, rigid plastic currettes for uterine aspiration. Moulded markings for size and correct tip orientation. Seven sizes available in either curved or straight designs. Connects to Pennine rotary swivel handle and wide bore, high vacuum tubing. Available as sets pre-attached to 2.0m thick walled aspiration tubing.

Aspiration Currettes are also available in boxes of 50, simply add '/50' suffix to product code.

Product Code	Type	Qty. Per Box	Product Code	Qty. Per Box
Curved Tip			2m Sets	
KC-1007	7mm	10	KCS-2007	10
KC-1008	8mm	10	KCS-2008	10
KC-1009	9mm	10	KCS-2019	10
KC-1010	10mm	10	KCS-2010	10
KC-1011	11mm	10	KCS-2011	10
KC-1012	12mm	10	KCS-2012	10
KC-1014	14mm	10	KCS-2014	10
Straight Tip			2m Sets	
KC-1507	7mm	10	KCS-2507	10
KC-1508	8mm	10	KCS-2508	10
KC-1509	9mm	10	KCS-2509	10
KC-1510	10mm	10	KCS-2510	10
KC-1511	11mm	10	KCS-2511	10
KC-1512	12mm	10	KCS-2512	10
KC-1514	14mm	10	KCS-2514	10

Karman Currettes

Clear, flexible currettes for uterine aspiration. Available with connector to fit to Pennine rotary swivel handle.

Product Code	Diam. mm	Qty. Per Box
AC-3004	4mm	50
AC-3005	5mm	50
AC-3006	6mm	50
AC-3007	7mm	50
AC-3008	8mm	50
AC-3009	9mm	50
AC-3010	10mm	50
AC-3011	11mm	50
AC-3012	12mm	50


Wide Bore Aspiration Tubing (Flexi-Cuff)

- Sterile, double wrapped individual peel pouches
- Leak proof seals
- Easily connects to currettes and Karman catheters
- Loose male connector in each pack for connection to suction machines where required
- Made of high density polythene
- Flexi-Cuff provides extra flexibility at a crucial point offering ease of use and improved technique

Product Code	Type	Qty. Per Box
KHS-3000	2m wide bore + rotary handle	10
KHS-3030	3m wide bore + rotary handle	5

Aspiration tubing is also available without the swivel handle:

Product Code	Type	Qty. Per Box
AT-2000	2m wide bore	10
AT-3000	3m wide bore	5

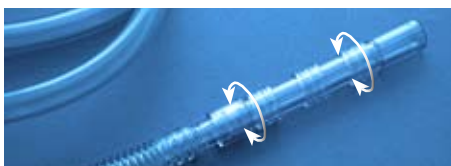

Please note: All angles and lengths are nominal

Tel: 44 (0) 1332 794880 Email: marketing@penninehealthcare.co.uk www.penninehealthcare.co.uk




Hy-G-Vac Suction Collection System

A complete range of compatible Pennine products designed for collection of surgical suction waste during suction procedures.


The core of the suction collection system is the Pennine Hy-G-Vac collection bottle:

- A safe, non-shatter, polypropylene bottle
- Non-removable leak-proof lid
- Quick and easy to assemble and use
- One product for all applications
- Hydrophobic block filter in the lid
- 2 litre capacity
- Safe transport after use including safe incineration properties
- HY-G-Sorb gelling agent to solidify bottle contents before disposal
- Range of trolleys, brackets and other devices to adapt and fix to existing hospital equipment

Product Code	Description	Qty. Per Box
HGV-2000	Hy-G-Vac suction collection bottle 2 litre capacity	20
HGV-2010	Link (tandem) connecting tube, F/F, 30cm, 6mm	100
HGV-2044	4-up floor trolley for HGV-2000	N/A
HGV-2042	2-up floor trolley for HGV-2000	N/A
HGV-2020	Specimen collection adaptor	10
HGV-2040	Ring holder for bracket or rail attachment	1
HGV-2055	Hy-G-Sorb gelling agent, with funnel	20
HGV-2050	HY-G-SORB gelling agent (20gm to gel approx. 1000ml water)	1
HGV-2052	HY-G-SORB gelling agent (4.0kg container)	1

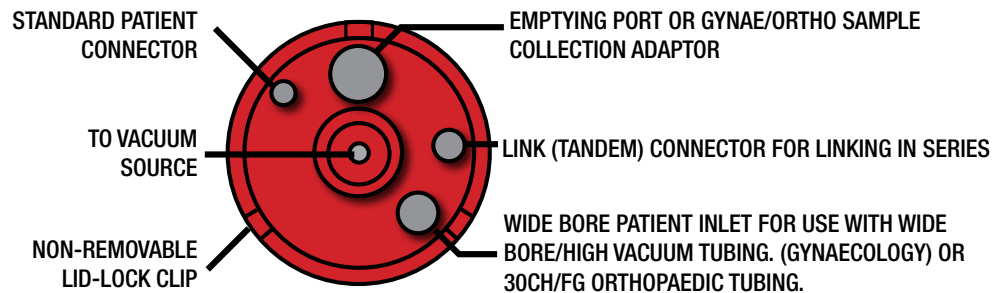

2-WAY, 3-WAY OR 4-WAY LINK UP SYSTEM (LINK PORT TO PATIENT PORT)

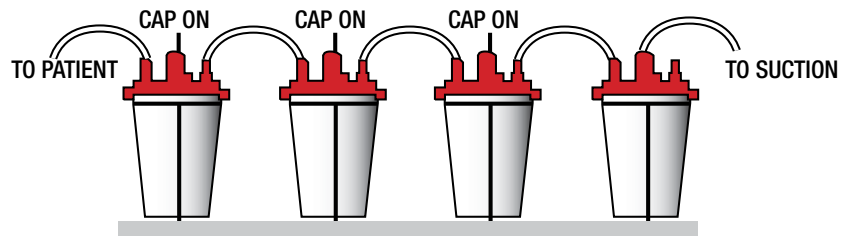



- Sterile and disposable
- Adhesive sharps disposal pad
- Sealable

Sharps Pad

For the containment of sharps prior to disposal. Enables used sharps to be counted, controlled and disposed of safely

Product Code	Description	Length/Size	Qty. Per Box
AP-1000	Sharps Pad	120mm x 70mm	100

- Sterile and disposable
- Adhesive backing for attachment to drape or trolley cover

Tip Cleaner

For the removal and cleaning of cauterised tissue from a diathermy pencil.

Product Code	Description	Length/Size	Qty. Per Box
TC-1000	Tip Cleaner	45mm x 45mm	300
TC-1000/100	Tip Cleaner	45mm x 45mm	100

- Standard switch button
- 3m connecting cable
- 2.5cm removable blade
- Supplied in sterile peel pouches

Diathermy Pencil

Disposable, mono-polar electro-surgical pencil for soft tissue surgery.

Product Code	Description	Qty. Per Box
EP-0001	Diathermy Pencil	20
EP-0101	Set with pencil, tip cleaner and holster	10

- Moulded PVC ring
- Supplied sterile, double wrapped
- Two strips of adhesive tape

Camera Drape

Clear drape, to be pulled over a video camera/light cable during preparation of the operating theatre.

Product Code	Description	Length/Size	Qty. Per Box
TP-0808/S	Camera Drape	2.5m	20

Pennine Healthcare is one of the UK's largest privately owned manufacturers of single-use sterile medical devices, with a customer base, not only in the UK but in over 50 countries worldwide.

Pennine has an extensive range of disposables for use in urology, anaesthetics, respiratory, gastroenterology and surgical suction, as well as being a leading supplier of Custom Procedure packs.

Pennine's vision is to be a "world class, innovative manufacturer, supplying the global healthcare market with outstanding commitment to its people, customers and suppliers".

Pennine is enormously proud of its achievements but knows that it is only successful because of its people and real values, which are at the heart of everything they do

To be progressive
To be responsive
To be passionate
To have integrity

With a culture and history of innovation and engineering expertise, Pennine is still a competitive supplier of disposables on a worldwide basis and today embraces the skills and techniques of 'Lean Sigma' manufacture, ensuring that efficiency and waste reduction are targeted and that Pennine staff are highly trained, motivated and entirely focussed on the customer.

Pennine is consistently pursuing new opportunities and embracing innovative solutions, working with clinicians and healthcare professionals to develop new products for advanced patient care.

For further information, please call our sales team on +44 (0) 1332 794880 or email us at: marketing@penninehealthcare.co.uk. Alternatively, you can view all our product ranges at our website www.penninehealthcare.co.uk



Gastroenterology



Urology



Custom Procedure Packs



Respiratory Products



Microsuction Cannulae



Airways Management

Pennine[®]
H E A L T H C A R E

Pennine Healthcare

City Gate
London Road
Derby DE24 8WY
United Kingdom
Telephone: +44 (0) 1332 794880
Fax: +44 (0) 1332 794890

Email: info@penninehealthcare.co.uk
Web: www.penninehealthcare.co.uk

Some of the products in this catalogue are subject to a minimum order quantity.
Please ask our sales department for further details.

Ivor Shaw Ltd trading as Pennine Healthcare