

Pelican Healthcare Ltd
Cardiff Business Park
Cardiff
CF14 5WF

www.pelicanhealthcare.co.uk
e-mail: contactus@pelicanhealthcare.co.uk

ORDER LINE: 0800 013 1220

SAMPLE LINE: 0800 052 1394

pelican
FEMININE
HEALTHCARE

SUILIT/07/10

single-use
medical instruments





pelican
HEALTHCARE
an eakin company



PELijelly Lubricant

PELijelly is a water based lubricant designed for use in gynaecological, digital and instrument examinations.

The sachet format allows minimal waste and also removes any concerns regarding cross infection.

-  Sterile 5g sachet
-  No effect on DNA
-  Silk like consistency to minimise discomfort from procedure
-  100 pack dispensing carton
-  Convenient sachet presentation
-  Fragrance free

Code	Description	Pack Size
400400	PELijelly Lubricant	100

Pelican PELIpack

The Pelican PELIpack contains all instrumentation necessary to either remove and fit an intra-uterine device or simply just remove.


This procedure pack will allow ease of use and save time on the set up of individual devices to perform the procedure.

All instruments from the Pelican PELIpack are also available separately.



-  Choice of PELIpacks to minimise waste and ease of use
-  Choice of PELIspec size on both packs
-  Time saving complete procedure pack
-  Instruments mirror design of metal equivalents for ease of use
-  Sterile and single use
-  Unique product identification number for traceability


Code	Description	Pack Size
400351	PELIpack IUD removal & fitting kit Medium PELIspec	4
400352	PELIpack IUD removal & fitting kit Medium Long PELIspec	4
400361	PELIpack IUD removal kit Medium PELIspec	4
400362	PELIpack IUD removal kit Medium Long PELIspec	4

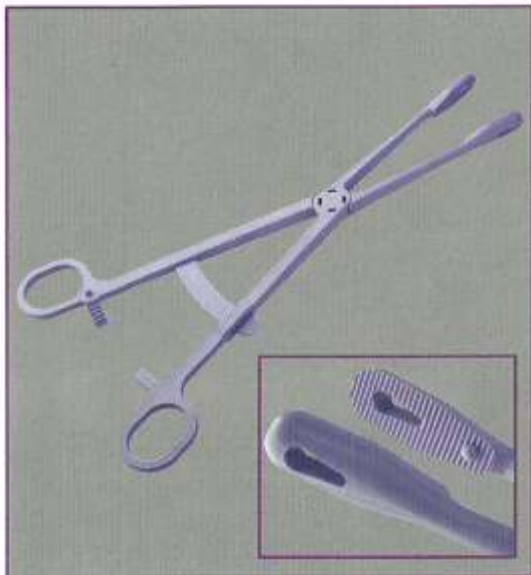



Pelican Spongeholder


The Pelican Spongeholder is designed to ensure visibility is not impeded when used in conjunction with the PELIspec.


The Spongeholder has been improved to demonstrate greater instrument strength and the secure grip of swabs.


 Secure tip grip for maximum hold

 Mirrors the metal equivalent for ease of use

 Enhanced secure locking mechanism

 Unique product identification number for traceability

 Specifically designed for use with swabs


 Choice of pack size


Code	Description	Pack Size
400300	Spongeholder	50
400302	Spongeholder	5


Pelican Dilator


The Pelican Dilator is designed to gently dilate the cervix, most commonly for the introduction of an intra-uterine device.


The tapered Pelican Dilator allows gradual dilatation from the two size options allowing 1mm increments from 2mm - 6mm.

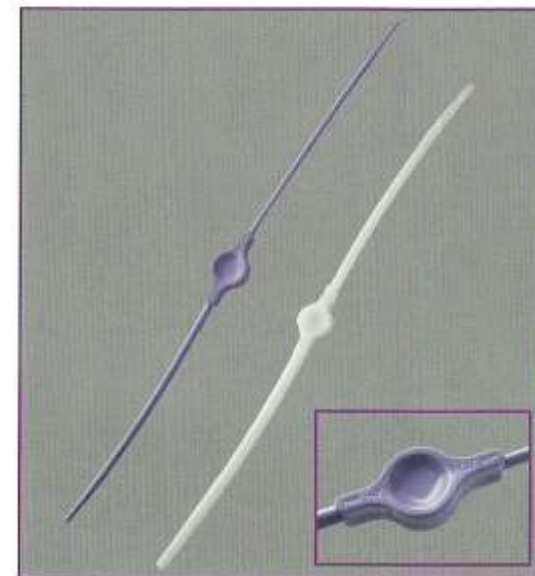



 Central thumb hold for ease of use

 Smooth rounded ends to minimise discomfort


 Unique product identification number for traceability

 Specifically designed for cervical dilation

 Tapered ends allow gradual cervical dilatation




 Available in two sizes

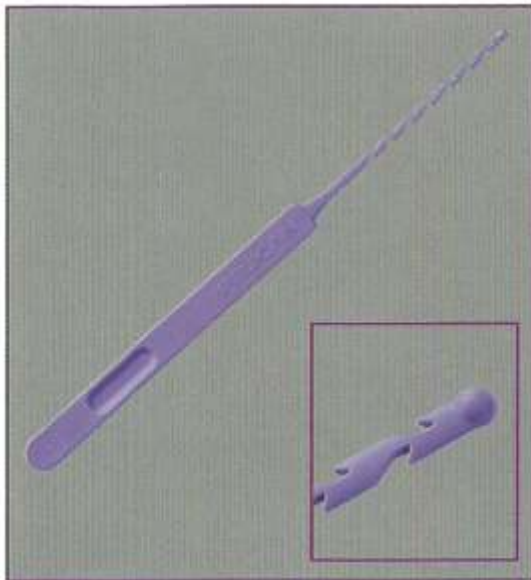Code	Description	Pack Size
400342	Dilator - 3/4mm	100
400343	Dilator - 3/4mm	10
400344	Dilator - 5/6mm	100
400345	Dilator - 5/6mm	10



Pelican Thread Retriever

The Pelican Thread Retriever is designed to locate lost IUD threads and can also be used to sound the uterus for the presence of the intra-uterine device.

-  Thumb recess for enhanced grip and control
-  Flexible stem to sound the uterus
-  Specifically designed notches for accurate location of the IUD threads


-  Seven notches on either side at 10mm intervals for optimum retrieval
-  Long handle (245mm) for ease of use and greater control
-  Smooth rounded tip to minimise discomfort

Code	Description	Pack Size
400560	Thread Retriever	5

Pelican Coil Remover

The Pelican Coil remover is specifically designed to remove intra-uterine devices and is characterised by an angled tip with a 3:2 tooth combination.

The tip of the device has a 'V' design where it is essential the coil threads are located (see inset picture) for a smooth retrieval. Specific instructions for use are included with every pack.

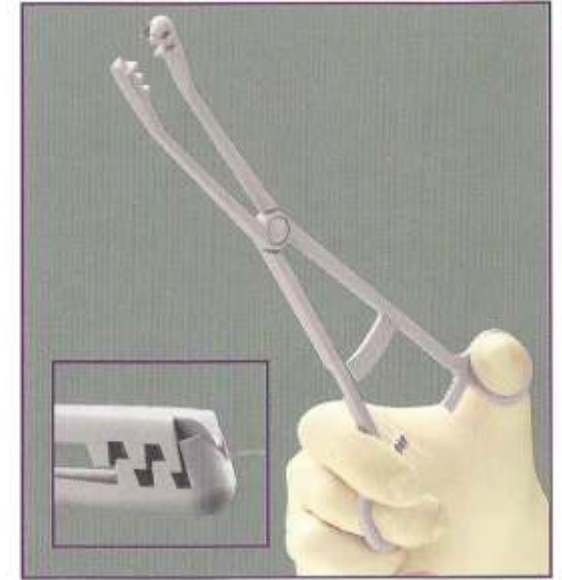

-  Secure lock mechanism to ensure a firm hold
-  Angled tip ensuring greater accuracy in locating coil threads
-  Carefully designed tooth combination to ensure accuracy of removal

-  Specifically designed instrument for ease of use
-  Designed for ease of use within the PELIspec without reducing visibility
-  Unique product identification number for traceability




Code	Description	Pack Size
400317	Coil Remover	50
400318	Coil Remover	5




Pelican Uterine Sound

The Pelican uterine sound is specifically shaped to aid insertion.

A slender flexible gauge with indicators which denote uterine depth prior to insertion of an intra-uterine device.

-  Notches allow measurement of the uterine depth
-  White in colour to aid visualising the depth
-  Mirrors the metal equivalent in design for ease of use


-  Increased patient comfort due to design
-  Unique product identification number for traceability
-  Manufactured in the UK

Code	Description	Pack Size
400310	Uterine Sound	100
400311	Uterine Sound	10




Pelican Scissors

The Pelican scissors are available in both a long and short handled version.

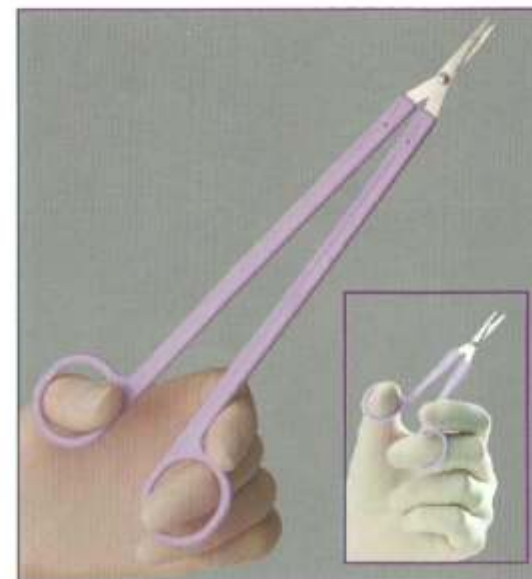
The long handled scissors have been specifically designed to use in conjunction with the PELIspec and give greater control and accuracy when cutting the coil threads.



-  Long handled option specifically designed to cut coil threads
-  Sharp blade
-  Firm finger grip for ease of use

-  Unique product identification number for traceability
-  Choice of scissor size
-  Choice of pack size

Code	Description	Pack Size
400320	Long Scissors	50
400321	Long Scissors	5
400322	Short Scissors	100

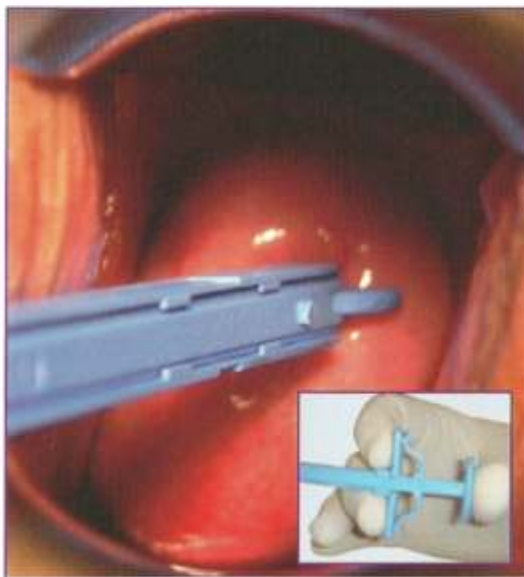


Cervix Biopsy Forceps

The Cervix Biopsy Forceps has been specially designed to obtain a biopsy sample with efficiency, safety and simplicity.


The Cervix Biopsy Forceps has a simple yet effective piston handle to ensure an adequate sample is obtained.

The ergonomic design allows the user full control to capture the sample tissue fragment.

 Innovative thumb pressure technique

 Sterile and single-use

 Ergonomic design for secure grip and maximum control

 Easy to use to obtain an adequate sample

 Latex free

 Effective piston sampling mechanism

Code	Description	Pack Size
01.390	Cervix Biopsy Forceps	25

Pelican Proctoscope

The Pelican Proctoscope is sterile and designed to minimise patient discomfort on insertion. The use of PELLjelly lubricant is recommended to reduce discomfort further during the procedure.

The Pelican Lighting system can be used in conjunction with the Pelican Proctoscope. The insertion point for the light source is sealed to ensure no risk of cross-infection during its usage.



 Available in 3 sizes

 Sterile


 Unique product identification number for traceability

 Crystal clear

 Latex free

 Manufactured in the UK

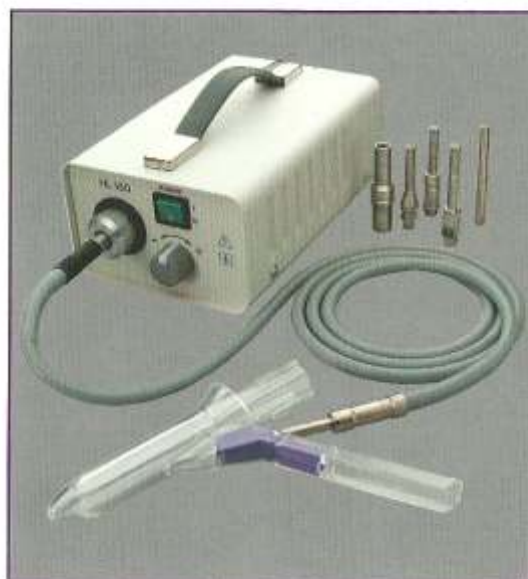


Code	Description	Pack Size
400370	Pelican Proctoscope - Paediatric	25
400380	Pelican Proctoscope - Medium	25
400390	Pelican Proctoscope - Large	25






Pelican Lighting System

A complete fibre optic lighting system is available to fit the Pelican Proctoscope range of instruments.

A range of adaptors are also available to allow connection of your current light source to the Pelican Proctoscope.

-  1.8 metre fibre optic illumination cable
-  150 watt Halogen light source
-  Detachable adaptor on cable to allow sterilisation

-  Output intensity adjustment
-  12 months warranty
-  Conforms to European standards

Code	Description	Pack Size
400500	Pelican Fibre Optic Lighting System	1
400510	Naughton Morgan/Storz Olympus Adaptor	1
400520	Heine Adaptor	1
400530	Woolf Adaptor	1
400540	Welch Allyn Adaptor	1
400550	Sapimed Adaptor and Cable	1

Description	Pack Size	Order Code
PELjelly Lubricant	100	400400
PELlpack IUD removal & fitting kit - Medium PELIspec	4	400351
PELlpack IUD removal & fitting kit - Medium Long PELIspec	4	400352
PELlpack IUD removal kit - Medium PELIspec	4	400361
PELlpack IUD removal kit - Medium Long PELIspec	4	400362
Spongeholder	50	400300
Spongeholder	5	400302
Dilator - 3/4mm	100	400342
Dilator - 3/4mm	10	400343
Dilator - 5/6mm	100	400344
Dilator - 5/6mm	10	400345
Vulsellum	50	400335
Vulsellum	5	400336
Tenaculum	50	400330
Tenaculum	5	400331
Thread Retriever	5	400560
Coil Remover	50	400317
Coil Remover	5	400318
Uterine Sound	100	400310
Uterine Sound	10	400311
Long Scissors	50	400320
Long Scissors	5	400321
Short Scissors	100	400322
Cervix Biopsy Forceps	25	01.390
Proctoscope Paediatric	25	400370
Proctoscope Medium	25	400380
Proctoscope Large	25	400390
Pelican Fibre Optic Lighting System	1	400500
Naughton Morgan/Storz Olympus Adaptor	1	400510
Heine Adaptor	1	400520
Woolf Adaptor	1	400530
Welch Allyn Adaptor	1	400540
Sapimed Adaptor and Cable	1	400550